


Informational Bulletin

IMPORTANT INFORMATION – PLEASE READ AND KEEP

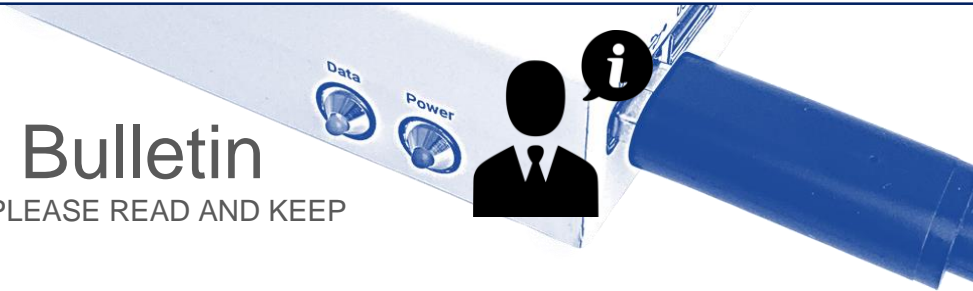

ARKON Connection to iMDsoft - MetaVision

- Document:** 063_ARKON_iMDsoft_MetaVision_IB-03_A-00_EN
- Date:** 2020-04-08
- Author:** Stephan Lang
- Related version:** All
- Distribute to:**
- Internal Use Only (Confidential)
 - Medcaptain
 - Dealer
 - End-user



Connection with ARKON possible

For the connection of the ARKON to the iMDsoft MetaVision PDMS **no new ARKON variant is required**. The connection is possible with the ARKON.

The data of the MIRUS - MACset, etCO₂, etVA as well as the alarms, are transmitted via ARKON to the patient monitor. Patient monitors in turn are connected to the PDMS, i.e. also to iMDsoft MetaVision, and transmit all available data. Consequently, also the transmitted data of the MIRUS.



Therefore, we do not need a direct connection of the ARKON to the iMDsoft MetaVision, but transmit directly to the patient monitor. **At the moment this is limited to the Philips IntelliVue, i.e. the use of the ARKON-PH.** Also with future ARKON variants the PDMS connection of iMDsoft MetaVision be possible in this way.